
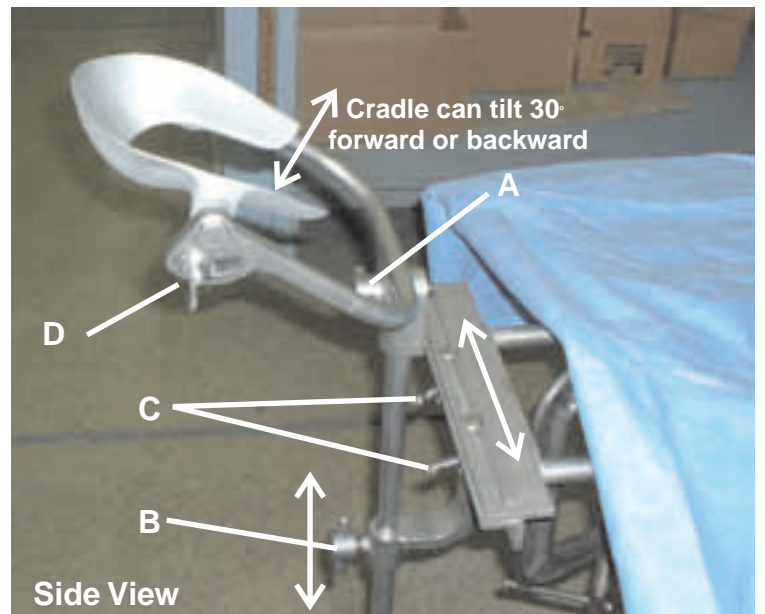
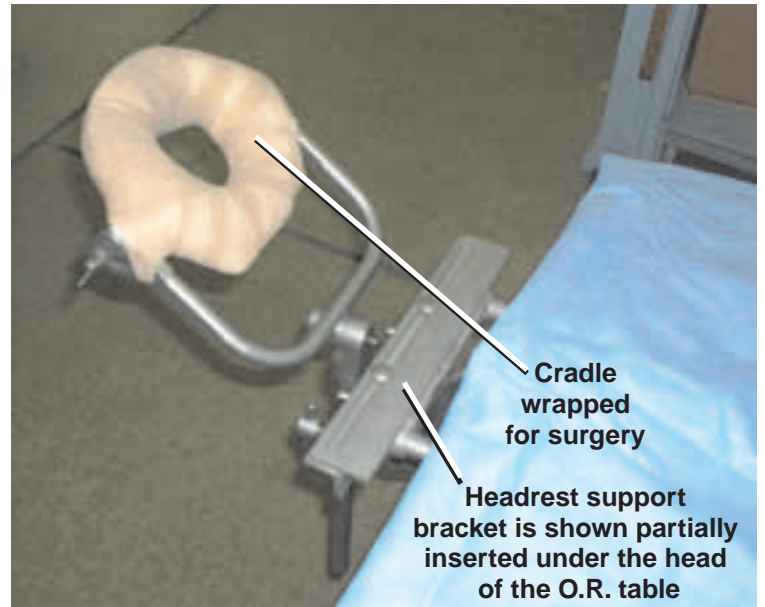


MODIFIED KING'S HEADREST

Indications:

The KAPP King's Headrest is used to position the patient's head during all surgical procedures. During surgery, the headrest can position and stabilize the patient's head at various angles relative to the configuration of the surgical table.

- **Most Durable**
- **Most Versatile**

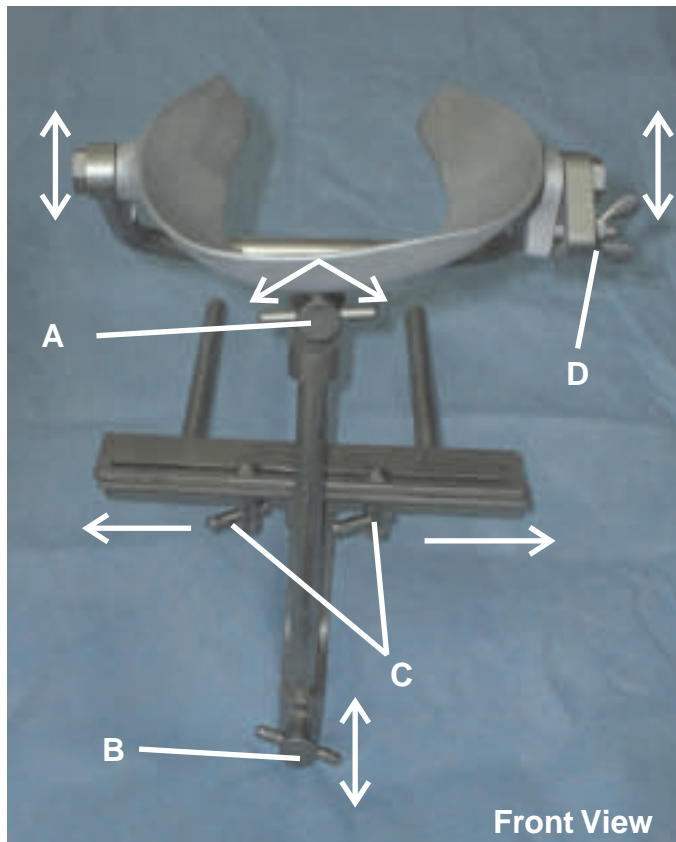

Features:

- **Thumb screw adjustment points allow a wide range of fixed positioning and stabilization:**
 - A - Adjust patient's head to turn 180° side to side.
 - B - Adjust headrest to move up and down (ceiling to floor) approx. 15 inches.
 - C - Adjust headrest to move medially or laterally (3.5" each way) to center line of table.
 - D - Adjust cradle of headrest to tilt patient's chin up or down approx. 30° in each direction.
- **The KAPP Modified King's Headrest fits all surgical O.R. tables.**
- **Manufactured from Surgical Grade 400 Series Stainless Steel.**
- **Lifetime Guarantee under normal O.R. use.**



KAPP SURGICAL INSTRUMENT, INC.

4919 Warrensville Center Road., Cleveland, Ohio 44128
(216)587-4400 FAX (216) 587-0411 1-800-282-KAPP
Email: info@kappsurgical.com Web Site: www.kappsurgical.com